

# Advance Advance Plus

Intermittent Catheters

## Designed to offer choices that empower life.

Developed with life's details in mind, Advance and Advance Plus intermittent catheters provide the choices that help people manage continence with confidence. Rest assured that whichever Advance catheter you choose, it's always the safe way forward.

### Safe Preparation

The unique gel reservoir ensures no additional lubrication is needed. Its self-contained gel supply reduces the inconvenience and risk of contamination from using an external source of water or other lubricant during preparation.


### Safe Handling

The 'no-touch' system helps avoid any transfer of bacteria from the surroundings to the catheter by preventing any direct contact with it. Instead, the catheter stays protected inside the sterile packaging up until the moment of insertion.


### Safe Application

The protective tip reduces the possibility of transferring bacteria from the distal region of the urethra further into the bladder. Bypassing this area helps reduce the risk of urinary tract infection.


### Protective Tip – for safe insertion

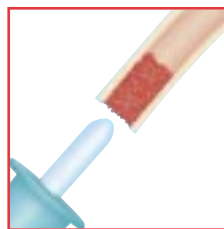

In most people, the greatest concentration of bacteria is located within the first 15mm of the distal urethra.

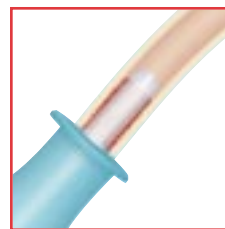

As the unique protective tip is inserted into the urethra, the tip shields the sterile catheter from contact with these bacteria.

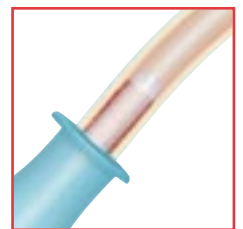


When forwarded into the urethra, the sterile catheter bypasses the concentration of bacteria, reducing the risk of pushing the bacteria further into the urinary tract.


Laboratory tests comparing standard hydrophilic catheters with the Advance intermittent catheter were run to measure bacterial transfer onto the catheter during handling.

The experiments tested for *Staphylococcus aureus* and *Escherichia coli*. The results clearly showed that the Advance intermittent catheter, which uses the "no-touch" system, had significantly less contamination.

Relative Recovery of *Staphylococcus aureus*

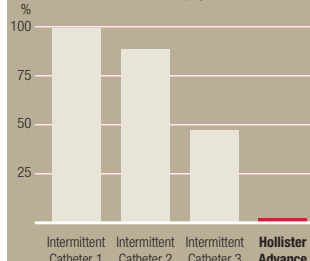

Relative Recovery of *Escherichia coli*

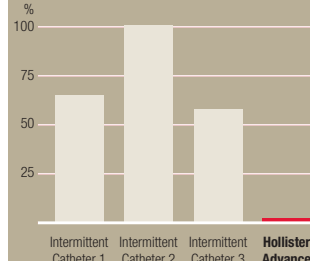


# Advance

## Intermittent Catheter

Selecting the Advance intermittent catheter helps bring more safety and convenience to the people you help.



- 1. Protective Ring Cap**  
for safety and ease of use
- 2. Unique Gel Reservoir**  
ready to use without water or additional lubrication
- 3. Unique Protective Tip**  
helps reduce the risk of urinary tract infection
- 4. Ultra-smooth Eyelets**  
help minimise potential trauma to the urethra

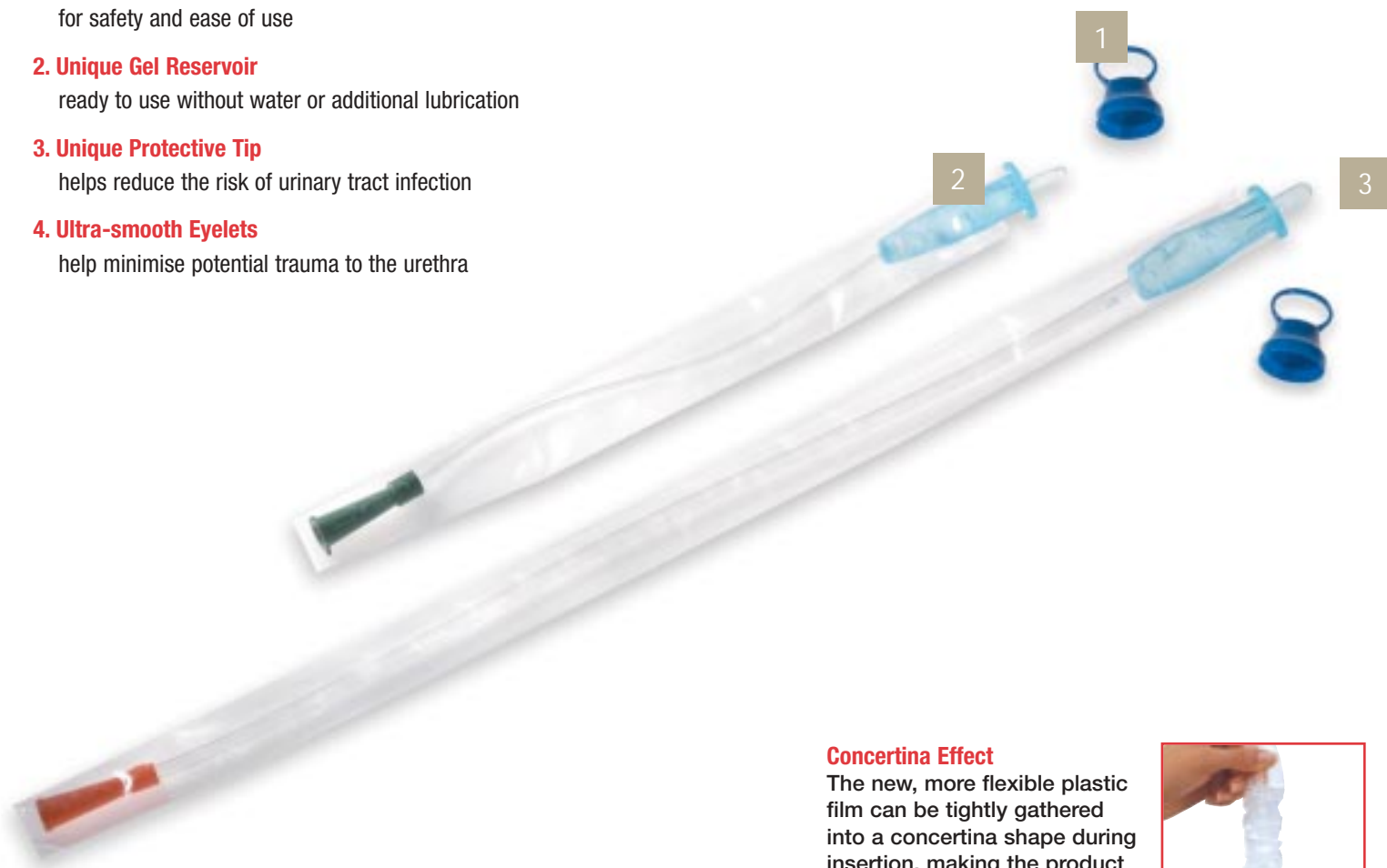

### Concertina Effect

The new, more flexible plastic film can be tightly gathered into a concertina shape during insertion, making the product easy to manipulate, even with one hand.



### Simplified Opening

Both products are supplied with finger holes in the packaging to provide easier opening, in a way that even those with limited manual dexterity will find simple.



### Simplified Hygiene

The protective ring cap helps maintain the sterility of the protective tip once the product is opened. It is designed to be easy to remove and replace.

# Advance Plus

## Intermittent Catheter

For people seeking independence, the Advance Plus intermittent catheter provides a safe option that can be used anytime, anywhere.




### Mobility Plus


The integrated 1500ml bag allows discreet and secure urine collection, whether indoors or out and about. The unique gel reservoir, together with the ring cap, prevents leakage by sealing off the bag.

### Handling Plus

The built-in handle makes the product easier to hold and to empty.

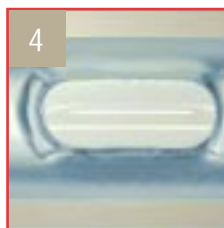

Ultra-smooth rounded eyelets formed with our advanced technology

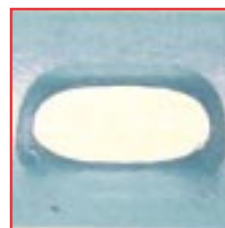

Close-up of eyelets produced on other catheters using traditional processes

### Simplified Insertion

The unique gel reservoir ensures that the catheter is evenly lubricated as it travels through the reservoir. Users can also dispense extra gel to meet their individual needs simply by squeezing the reservoir during insertion. The ability to regulate the lubrication, together with the water-soluble gel, helps provide maximum comfort during insertion and removal of the catheter.

### Simplified Protection

The catheter's distinctive eyelets are ultra-smooth, improving comfort during insertion and removal. Their rounded shape helps protect sensitive urethral tissue.

# A relationship that makes a difference

## **Building Relationships**

As a healthcare professional, you develop a trusted bond with the people you help, returning them to their families and their lives. You also link Hollister to the lives of these people. Through listening to you and to them we develop solutions that can help others every day.

## **Sharing Goals**

As an expert in helping people manage their day-to-day health, you strive to do the very best. And Hollister aims to help by providing more than just the right products. With us you'll find a range of support materials and services available. Our experienced Customer Service Team have been highly trained to provide expertise and advice, and our educational support materials assist you when training people how to catheterise. All to help you achieve your goals, and to feel confident those you help can also achieve theirs.

## **Focusing on Details**

Sometimes the smallest things can also be the most important. At Hollister we focus on the details because we see them as vital to building the big picture. And it's something you'll see reflected in all our products, with features like the anti-reflux flap on our urinary sheaths that prevents backflow and increases wear time, or the unique catheter protective tip that helps reduce the risk of urinary tract infections. By concerning ourselves with these details, and by carefully listening and paying attention to people's needs at all levels, we can provide every possible advantage for you and those you help.

## **Delivering Confidence**

For more than 80 years, Hollister has been making a difference where it matters, playing a role in making life rewarding and dignified. As we continue to develop products and services, our goal remains the same – to promote independence and quality of life.



<b>Advance</b> Intermittent Catheter	Size	Length	Tip	Colour Code	Stock No.	Box Quantity
	Ch06	20cm	Nelaton	●	92062	25
	Ch08	20cm	Nelaton	●	92082	25
	Ch10	20cm	Nelaton	●	92102	25
	Ch12	20cm	Nelaton	○	92122	25
	Ch14	20cm	Nelaton	●	92142	25
	Ch08	40cm	Nelaton	●	92084	25
	Ch10	40cm	Nelaton	●	92104	25
	Ch12	40cm	Nelaton	○	92124	25
	Ch14	40cm	Nelaton	●	92144	25
	Ch16	40cm	Nelaton	●	92164	25
	Ch18	40cm	Nelaton	●	92184	25
	<b>Advance Plus</b> Intermittent Catheter	Size	Length	Tip	Colour Code	Stock No.
	Ch06	40cm	Nelaton	●	94064	25
	Ch08	40cm	Nelaton	●	94084	25
	Ch10	40cm	Nelaton	●	94104	25
	Ch12	40cm	Nelaton	○	94124	25
	Ch14	40cm	Nelaton	●	94144	25
	Ch16	40cm	Nelaton	●	94164	25
	Ch18	40cm	Nelaton	●	94184	25
	Ch12	40cm	Tiemann	○	95124	25
	Ch14	40cm	Tiemann	●	95144	25
	Ch16	40cm	Tiemann	●	95164	25

Advance and Advance Plus catheters are available on Drug Tariff.

Advance and Advance Plus intermittent catheters are part of a range of Hollister products designed specifically to help people manage their continence.

If you'd like to receive a sample,  
please call our **FREEPHONE NUMBER 0800 592202**.  
Customers in Ireland please call **FREEPHONE NUMBER 1800 503400**.

**Hollister Limited**

Rectory Court, 42 Broad Street, Wokingham, Berkshire RG40 1AB, England.  
Freephone: 0800 521377 Tel: 0118 989 5000 Fax: 0118 989 5060

**Hollister Limited**

Belgard Road, Tallaght, Dublin 24, Ireland.  
Freephone: 1800 503400 Tel: 01 404 1680 Fax: 01 404 1682

[www.hollister.co.uk](http://www.hollister.co.uk)

# Advance Advance Plus

Intermittent Catheters



Choices that empower life

