

The right choice.

Ordering Information

Stock No.	Description	Quantity
9799	Anchor Fast Oral Endotracheal Tube Fastener	Box of 12
9788	Anchor Fast Replacement Neck Band Straps	Box of 4

References

- Veneman W, Adams-Stubbs D. (2010) An Alternative Approach to Prevent Lip Pressure Ulcers in the Intubated Patient. Presented at The Symposium on Advanced Wound Care. 17-20 April 2010, Orlando, Florida.
Wade Veneman, RCP, RRT and Debbie Adams-Stubbs, RCP, RRT. Community Regional Medical Center, Fresno, CA USA. (2008-2010) Number of subjects: unknown, observational study of all ventilated patients from November 2008 to February 2010.
- Miller K, Fenstermaker C, Limoge D, Leshko R, Lutz A, Morgan G, Cornman L. (2009) A Utilization of Anchor Fast Oral Endotracheal Tube Fastener to Reduce the Incidence of Lip Ulcers. Presented at The 55th International Respiratory Congress. 5-8 December 2009, San Antonio, Texas.
Kenneth Miller MEd, RRT-NPS, Chris Fenstermaker, BS, RRT-NPS, Diane Limoge CCRN, Robert Leshko BS, RRT, Angela Lutz, BS, RRT-NPS, Georgiann Morgan, Linda Cornman, BS, RRT-NPS. Lehigh Valley Health Network, Bethlehem, PA USA (2008-2009) Number of subjects: unknown, observational study of all ventilated patients from January 2008 to December 2008.
- Garcia R, Jendresky L, Colbert L, Bailey A, Zaman M, and Majumder M. Reducing Ventilator-Associated Pneumonia Through Advanced Oral-Dental Care: A 48-Month Study. American Journal of Critical-Care. 2009. <http://www.ajconline.org>



Secure
Clean
Practical



Hollister and logo, and Anchor Fast are trademarks of Hollister Incorporated.
© 2011 Hollister Incorporated. C00208.

Hollister Limited
Rectory Court
42 Broad Street
Wokingham
Berkshire RG40 1AB

**For more information and
to request FREE samples:**
Tel: 0800 521 377
Email: samples.uk@hollister.com

www.hollister.com/uk


The right choice for securing endotracheal tubes.
The right choice for optimising oral care.
The right choice for helping safeguard skin health.



Anchor Fast

Oral Endotracheal Tube Fastener

The Anchor Fast Oral Endotracheal Tube Fastener combines a number of unique features that make it the right choice. As part of a multi-faceted approach, it helps to prevent the formation of lip ulcers^{1,2} and eases access to the oral cavity to optimise patient oral care. Implementation of a comprehensive oral and dental care protocol has been shown to help reduce the incidence of ventilator-associated pneumonia (VAP) during endotracheal intubation.³

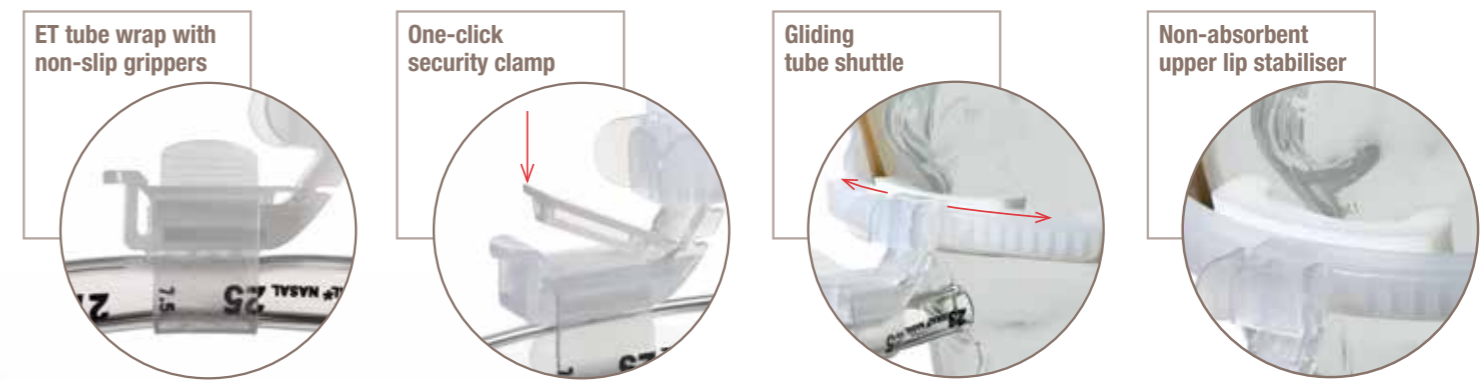

Simple to use

1. place

place non-absorbent upper lip stabiliser so it rests on skin between nose and lip; place skin barrier pads on cheeks

2. secure

secure using ET tube wrap with non-slip grippers and one-click security clamp

3. position

position gliding tube shuttle; reposition for oral care and pressure relief

One-click security clamp

locks the ET tube in place; one-click releases the clamp for tube repositioning

Adjustable neck strap

helps provide extra security against possible extubation and it is replaceable

Skin barrier pads

fasten the device firmly in place; the unique skin barrier formulation is designed to be gentle on skin, easy to apply and remove, and moisture resistant to help prevent the device from moving out of position

Non-absorbent upper lip stabiliser

keeps the tube lifted off of the patient's lip, which helps alleviate undo pressure on the patient's mouth and helps to prevent excess moisture build-up that can lead to skin breakdown

Gliding tube shuttle

allows the ET tube to be easily repositioned from side to side, while the tube remains secure

ET tube wrap with non-slip grippers

helps prevent inadvertant tube movement



Reducing VAP Through Advanced Oral-Dental Care

In an independent study, "findings suggest that use of advanced tools, a comprehensive oral care protocol, and staff compliance with the protocol can significantly reduce rates of ventilator-associated pneumonia and associated costs."*³

*Endotracheal tube fasteners were not a subject of the study.